

Tracking Chart

To better help you track your skin's progress, keep a record of when and how you used your DPL™ System and how your skin looked and felt afterwards.

10 Week DPL™ Treatment Plan

Week	# of Treatments	# of Minutes	Visible Results



Quick Start Guide

LED Technologies, LLC uses a proprietary array of Light Emitting Diodes or LEDs to produce a unique combination of frequencies and wavelengths which penetrates deep into the skin layers to help remove the signs of aging; fine lines, wrinkles, crow's feet and blotchy skin. The DPL™ Therapy System is extraordinarily effective at making your skin look younger and more healthy.

Start enjoying the benefits of "Deep Penetrating Light" Today!




1. Unpack your DPL™ System.

Your box will include the following:

- ◆ DPL™ System (2 panels)
- ◆ DPL™ Stand (3 parts)
- ◆ Wall Adapter
- ◆ Panel Connector Wire
- ◆ Quick Start Guide
- ◆ User Guide
- ◆ Warranty Card


2. Assemble Stand (3 pieces).

Insert column (1) into base (2) (one or two clicks).
 Insert mounting bracket (3) onto column (1) and slide to the right until it fastens (loud click).

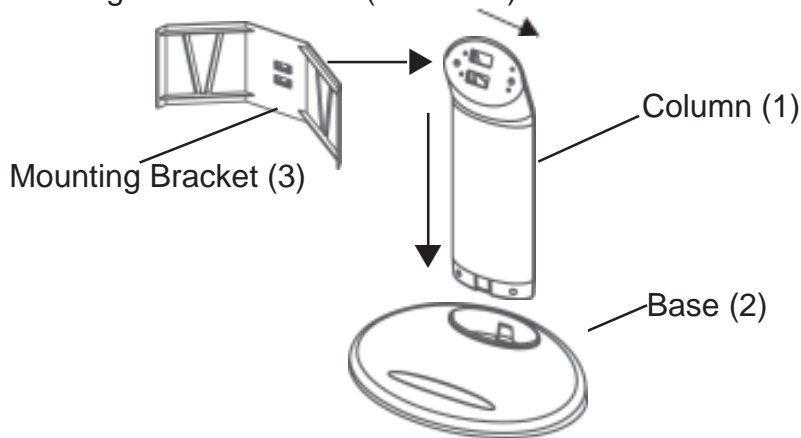


3. Attach DPL™ panels to stand.

The panels are different, the “Master” panel (A) has 2 plugs and a switch along the bottom edge.
 Slide Master panel (A) onto right side of mounting bracket.
 Slide Left panel (B) onto left side of mounting bracket.

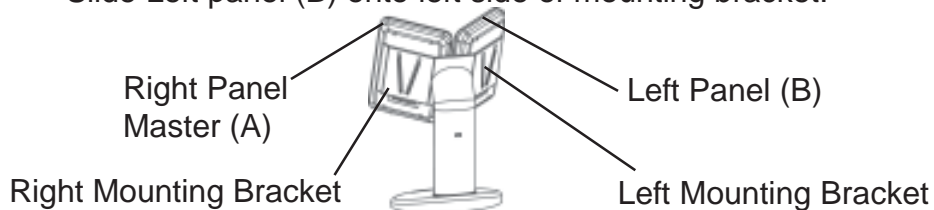


4. Attach electrical panel connector between panel A and panel B.

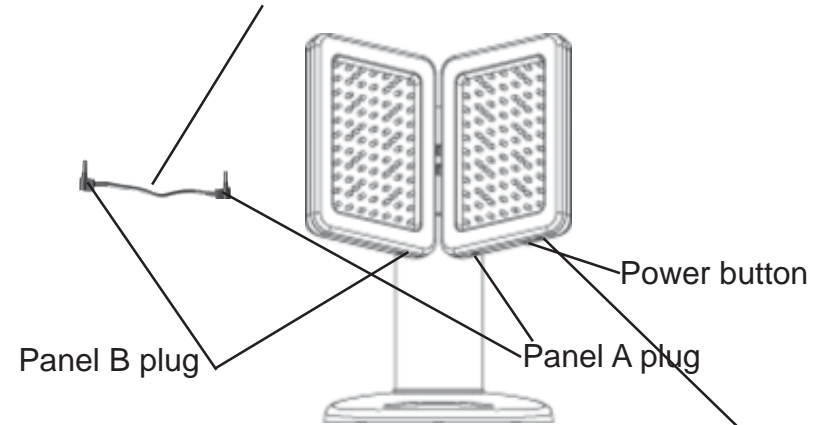

5. Attach power plug from power supply into Panel A.

6. Attach outlet plug into electrical outlet.

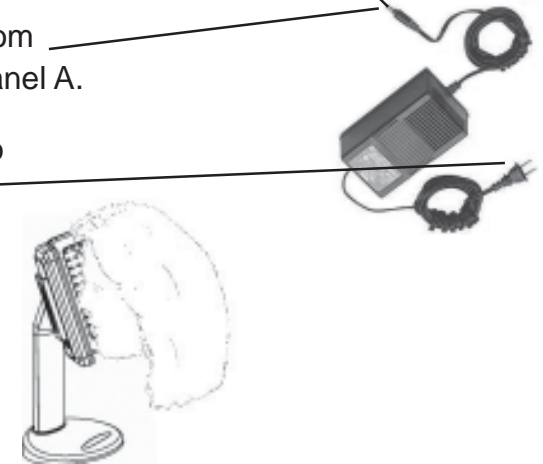

7. Place face approximately 1/2 inch from center of unit.

8. Press power button and release.

- ◆ You will hear a beep and the lights will come on.
- ◆ You can only see the Red LEDs; the others are in the infrared spectrum and can't be seen with the human eye.
- ◆ The unit will beep in 1 minute intervals for a total of nine minutes.

9. The unit will automatically turn off at the end of nine minutes.

Usage: The optimal amount of exposure is two (9 minute) treatments a day. Additional treatments will not improve your results.